




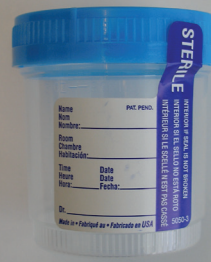

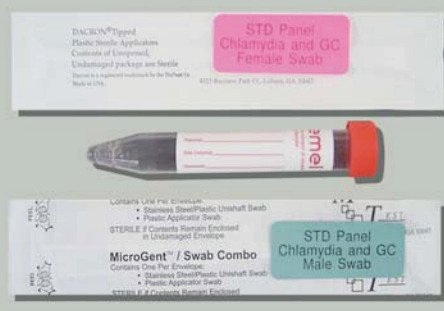




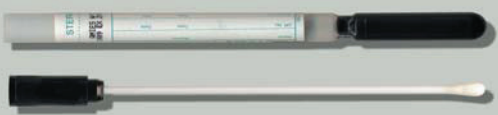


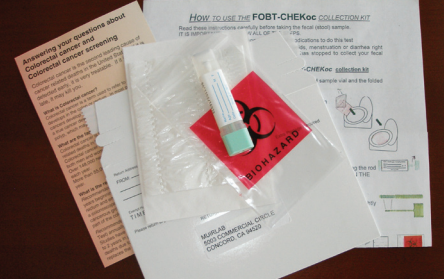





Dark red area represents Storage condition after collection.

TRANSPORT CONTAINER	TEST(S)	TRANSPORT CONTAINER	TEST(S)
	<b>Adult Blood Culture</b> 6 years + Aerobic Plus Order # 73067 Anaerobic Order # 48681 <b>Room Temperature</b>		<b>UA Collection Kit</b> • With Preservatives Order # 53450 <b>Room Temperature</b>
	<b>Pediatric Blood Culture</b> < 6 years Order # 1244 <b>Room Temperature</b>		<b>Urine C &amp; S Transport</b> Order # 62495 <b>Room Temperature</b>
	<b>Blood Culture for AFB or Fungus</b> Order # 47927 <b>Room Temperature</b>		<b>Sterile Cup</b> • Sputum - Transport immediately to Lab • Nasal Wash or Aspirate • Urine • Stool for H. Pylori Ag Order # 107135 <b>Sputum: Room Temperature</b> <b>Nasal, Urine, Stool: Refrigerate</b>
 Stuart's	<b>Routine Bacterial Culture Swab</b> • Aerobic Cultures • Rapid Strep Screen Order # 51767 <b>Room Temperature</b>		<b>STD Panel by PCR</b> • Chlamydia • GC Blue for male Order # 37419 Pink for female Order # 37418 <b>Refrigerate</b>
 Amies without charcoal	<b>Anaerobic Culture Swab</b> • Deep Wounds • Closed Abscesses Order # 51769 <b>Room Temperature</b>		<b>Affirm VPIII</b> • Bacterial Vaginosis/ Vaginitis Panel Order # 75506 <b>Room Temperature</b>
 Universal Viral Transport Swab Applicators	<b>UVT</b> (not for bacterial cultures) • Rapid Flu A & B • Chlamydia Culture • Rapid RSV • Mycoplasma Culture • Viral Culture • Ureaplasma Culture • HSV Culture Only • STD by PCR Order # 107583 <b>Refrigerate</b>		<b>For Stools</b> <b>C &amp; S Medium/Cary Blair</b> • Routine Culture • Clostridium Difficile • Giardia • Shiga Toxin (Cap colors may change) Order # 65823 <b>Refrigerate</b>
 Amies with charcoal	<b>Charcoal Culture Swab</b> • GC • Genital Cultures • Beta Strep Culture Order # 51768 <b>Room Temperature</b>		<b>Formalin &amp; PVA for Stools</b> • Ova & Parasite • Cryptosporidium • Cyclospora/Isospora Order # 66861 (Cap colors may change) <b>Room Temperature</b>
 Amies with charcoal with N/P Swab	<b>Charcoal Nasal Swab</b> • Pertussis Culture Order # 51771 <b>Room Temperature</b>		<b>Stool Collection Device Immunoassay</b> • Occult Blood <b>COLLECTION KIT</b> Order # 71031 <b>Room Temperature - 10 days</b>
 Nasopharyngeal Swab	• Pertussis PCR Order # 73337 <b>Refrigerated</b>		<b>ThinPrep Pap Test</b> • GC/Chlamydia Test with Pap Test • HPV Test with Pap Test Order #49865 TP Spatula & Brush Order #49864 TP Brooms Order #49863 TP Vials <b>Room Temperature</b>
 <b>MUIRLAB</b> PT Name: _____ DOB: _____ Collect Date ___/___/___ Time _____ Collector Name: _____	<b>Custom Label</b> • For patient safety, all specimens must be labeled. Order # 65805		